






The Porter AVS 5000 does **not retrofit** other manufacturers flowmeters. The PORTER Scavenger breathing circuit **does retrofit** the flowmeters listed below. Use these part numbers to *upgrade* the other manufacturer's flowmeters to the market leader.

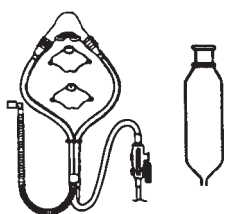


PORTER Retrofit Kits come complete with the following: #5155-3 (Adult) or #5155-4 (Pedo) Rubber Goods, #5501-RK Vacuum Hook-up Kit, #4100 3Liter Bag, and where applicable, any necessary bag tee or flowmeter adapters.

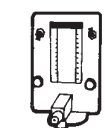

Accutron 5800-1A*

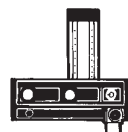

MDT-McKesson Analor Four and Compalor 5800-1A*

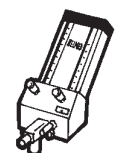

Belmed 5800-1A*



NCG Chemetron 5800-2A*




Coastal 5800-2A*




Parkell Sedetron 5800-1A*



Frazer Sweetman, Frazer Harlake, Ohmeda and Matrix 5800-3A*

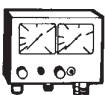


Spica 5800-7A*

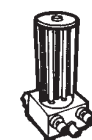


Hampton, Simplaire and Micro 5800-4A*



Starflite Denetron 5800-5A*

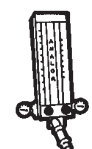

McKesson Analor I and II 5800-5A*

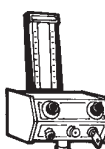


Trident 5800-1A*



McKesson Analor III and 3C 5800-6A*

PORTER Retrofit Scavenger

Replaces  Porter Clean-Air™ rubber goods **Non-Scavenging** (Mask with plastic sweeps) 5800-8A*

Replaces  Porter Clean-Air™ rubber goods and **Canister** Scavenging System 5800-9A*

*Pedo size inhaler substituted on request

IMPORTANT: It is not recommended to retrofit PORTER Scavenger rubber goods to Adec or Veriflo machines.

PORTER

NITROUS OXIDE SEDATION SYSTEMS



PORTER INSTRUMENT COMPANY, INC.


245 TOWNSHIP LINE RD. P.O. BOX 907 HATFIELD, PA 19440-0907 USA [215] 723-4000 / FAX [215] 723-2199
www.porterinstrument.com

FM-791 Rev. C 06/01



New AVS Technology... Exclusively From Porter

Our unique patent approved design automatically activates the scavenger vacuum as soon as gas begins to flow from the flowmeter. As a result, it prevents Nitrous Oxide/Oxygen from being administered without scavenging. In accordance with NIOSH guidelines, the vacuum should always be "ON" while administering Nitrous Oxide. The AVS 5000 retrofits to all flowmeters manufactured by Porter Instrument. The complete system now includes the flowmeter, the AVS 5000, and the scavenger breathing circuit.

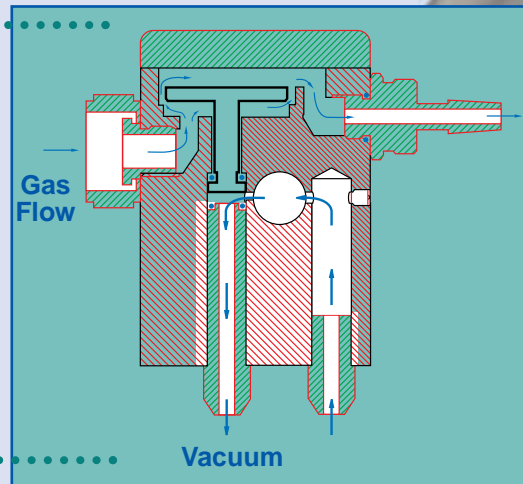

Convenient Operation –

The vacuum level sight glass lets the user monitor the levels of vacuum and make adjustments to the flow rate when necessary. Porter also designed the system to reduce tubing clutter.

When Gas Flows – Vacuum is "ON"... It's Automatic

When the procedure is finished, the vacuum is automatically shut down by the AVS. This is accomplished by turning off the flowmeter's Flow Control Knob.



Gray Autoclavable



Green Disposable*



White Citrus Disposable*

Worry-Free Scavenging System

- A flexible inner liner creates a tight seal around the nose for improved scavenging. Porter's unique double mask sets the industry standard for comfort and scavenging efficiency.
- Masks are autoclavable and available in Adult and Peds sizes. Disposable mask liners are also available for offices looking for added convenience.

Autoclavable or Disposable Liners –

**Disposable Liners are available in Unscented or Citrus scented.*



The Quality System for Porter Instrument Company, Inc. is Certified to ISO 9001. The scope of our registration is: "The design, manufacture, distribution and servicing of Dental Flowmeters, Gas Scavenging Systems, Dental Steam Sterilizers, Gas Distribution Systems and Office Communication Systems for use in the Dental Profession."